

## FOR THE ACTIVE WHEELCHAIR USER

### Trochanter support



firmer density to provide lateral stability | redistributes pressure to help create an ischial well

**Ischial well - Two compartments** for superior pressure reduction in the sacral and IT-coccyx area & extra positioning possibilities



**SmartCell – Foam barrier** limits sliding forward

**Ergonomically-shaped foam** for easy transfers

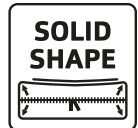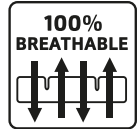

- ▲ Perfect for users who transfer frequently
- ▲ Excellent lateral and frontal stability
- ▲ Features a nifty storage pouch for personal items



The Vicair Active 02 is a pre-contoured cushion which consists of 4 individual compartments filled with Smartcells at the back of the cushion and an anatomically contoured foam inlay at the front. Striking a unique balance between preserving skin at high risk areas at the ischials, coccyx and sacrum by redistribution of pressures, stabilizing the pelvis and facilitating easy transfers. For specific needs the filling of each of the 4 SmartCell filled compartments can be easily adjusted. Allowing you to change the density of each compartment benefitting users who: (1) Need correction or accommodation of an asymmetric pelvis. (2) Want (partial) off-loading of the ischial, sacral and coccyx area. (3) Benefit from additional lateral support.

### AVAILABLE COVERS

Not applicable. The Vicair Active 02 only has one cover serving as inner and outer cover.

### PRODUCT PERFORMANCE

#### Skin protection / Anti decubitus

Pressure redistribution	▲▲▲▲△
Shear force reduction	▲▲▲▲△
Microclimate control	▲▲▲▲△
Impact & vibration damping	▲▲▲▲△

#### Positioning

Lateral stability	▲▲▲▲▲
Frontal stability	▲▲▲▲▲
Sliding reduction	▲▲▲▲▲
Allows manual adjustments	▲▲▲▲△

#### Maintenance & Safety

Durability	▲▲▲▲▲
Reliability	▲▲▲▲▲
Functionality when misused	▲▲△△△
Easy to maintain	▲▲▲▲△
Level of hygiene	▲▲▲▲△

#### Functionality & Application


Comfort	▲▲▲▲▲
Ultra lightweight	▲▲▲▲△
Allows easy transfers	▲▲▲▲▲
Available standard sizes	▲▲▲▲▲
Maximum user weight	▲▲▲▲△

### TECHNICAL SPECIFICATIONS

Thickness	6 and 9 cm
Average weight	1250 grams
Max. user weight	200 kg
Machine washable	60°C (without foam)

### CUSHION SIZES

- Standard (Sizes > 53 cm thickness 9cm only)
- On request

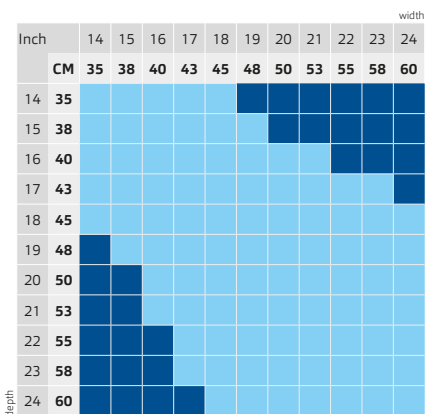

070069-OCT-EN-REV01/2025