

ICARE IC333 PRODUCT DETAIL

KEY FUNCTIONS




Head Lift



Knee Break



Height Adjustable
22-66cm

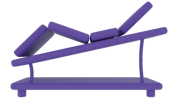

Trendelenburg and
Reverse Trendelenburg



BED CODES & SIZES

LONG SINGLE

Code: IC333LS

Mattress Surface Size: 2030 x 900

Base Size: 2150 x 960

KING SINGLE

Code: IC333KS

Mattress Surface Size: 2030 x 1070

Base Size: 2150 x 1120

LONG DOUBLE

Code: IC333D

Mattress Surface Size: 2030 x 1350

Base Size: 2150 x 1400

- One hand control

ONE PIECE QUEEN

Code: IC333Q

Mattress Surface Size: 2030 x 1520

Base Size: 2150 x 1580

- One hand control

SPLIT QUEEN (HEAD & KNEE FUNCTIONS) (ONE PIECE HI-LO & TRENDLENBERG)

Code: IC333SQ

Mattress Surface Size: 2030 x 750

Base Size: 2150 x 1580

- Individual hand control per side
- Fits a one piece queen headboard

SPLIT KING

For a split king order 2x IC333LS

Mattress Surface Size: 2030 x 900

Base Size: 2150 x 1920 (as partner)

KEY FEATURES

- 200kg Safe Working Load
- Optional Headboards and Footboards
- External Rechargeable Battery Back-up
- End and Corner Mattress Retainer Bars
- 10 Year Warranty

TECHNICAL SPECIFICATIONS

- Safe working load - 200kg on long single and king single sizes and 300kg for double and queen sizes.
- Hand Piece - Backlit for night time use. Easy press buttons. Wired stretch cord. Hook attachment for hanging to side of bed. Wipeable.
- Fabric - Hard wearing commercial grade fabric. Vinyl is treated with antibacterial properties.
- Mattress Surface - Strong metal mesh under mattress allowing increased airflow
- Actuators - 1x 6000 newton - 24 watt head adjust actuator. 1x 6000 newton - 24 watt knee adjust actuator. 2x 6000 newton - 24 watt height adjust actuators (2x 8000 newton on Double and Queen sizes)
- No pull motors for added safety.
- Wheels - 4x heavy duty swivel castors with two lockable castors.
- Head lift angle 62° and knee break angle 42°
- Max height - 66cm from floor to top of mattress platform.
- Min height - 22cm from floor to top of mattress platform.
- ARTG Identifier: 233386

COMPATIBLE ACCESSORIES AVAILABLE

- High Side Rail
- Low Side Rail
- Bed Stick
- U-Assist Rail
- Over Bed Pole
- Headboard / Footboard
- Slide-Out Side Rail
- Over Bed Tables
- Bedside Pouch
- Bedside Folding Shelf
- Bedside Tables
- Blanket Cradle