

icare

IC555

BARIATRIC

This is a bariatric model with a high weight rating. Featuring all functions and a 300kg safe working load.




KEY FUNCTIONS:



Head Lift



Knee Break

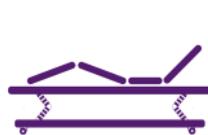

Height Adjustable
22-66cm



Trendelenburg and
Reverse Trendelenburg


COLOUR CHOICES:




Onyx



Stone

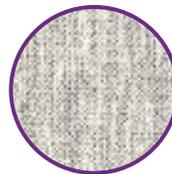

Silver

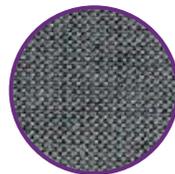

Smoke



Vinyl Ebony



Vinyl Dove

BED CODES & SIZES:



KING SINGLE

Code: IC555KS

Mattress Surface Size: 2030 x 1070

Base Size: 2150 x 1120

KEY FEATURES:



300kg Bariatric
Safe Working
Load



Strong Breathable
Mattress Platform



Ultra Low Design
Reducing Risk of
Falls



External
Rechargeable
Battery Backup



Lockable Heavy
Duty Castors



Available with
Headboards



End and Corner
Mattress Retainer
Bars



10 Year Warranty



Backlit Wipeable
Hand Piece



Unlimited Under
Bed Clearance

TECHNICAL SPECIFICATIONS:

- » Safe working load - 300kg safe working load. The safe working load includes any model of Icare mattresses.
- » Hand Piece - Backlit for night time use. Easy press buttons. Wired stretch cord. Hook attachment for hanging to side of bed. Wipeable.
- » Fabric - Hard wearing commercial grade fabric. Vinyl is treated with antibacterial properties and waterproof.
- » Mattress Surface - Strong metal mesh under mattress allowing increased airflow.
- » Actuators - 1x 6000 newton - 24 watt head adjust actuator. 1x 6000 newton - 24 watt knee adjust actuator. 2x 8000 newton - 24 watt height adjust actuators
- » No pull motors for added safety.
- » Wheels - 8x heavy duty swivel castors with one lockable castor on each corner.
- » Head lift angle 62°
- » Knee break angle 42°
- » Max height - 66cm from floor to top of base.
- » Min height - 22cm from floor to top of base.

ACCESSORIES AVAILABLE:

- » High Side Rail
- » Low Side Rail
- » Bed Stick
- » Over Bed Pole
- » Full Length Side Rails
- » Headboard
- » Footboard

