

Ocean VIP Ergo & Dual VIP Ergo

OUR MOST ADVANCED TILT-IN-SPACE SHOWER CHAIRS



ADVANCED PREMIUM PERSONAL CARE SOLUTIONS

CLINICALLY-LED DESIGN, ADJUSTABLE AND SUPPORTIVE



Yes, you can.®

AQUATEC®

Introducing our latest **Ocean** Dual VIP Ergo

Meet our most advanced tilt-in-space shower chairs with the latest design enhancements. Highly configurable and modular, the Ocean VIP Ergo and Ocean Dual VIP Ergo offer bathing and showering solutions with increased levels of comfort and support.



Ocean VIP Ergo

An advanced tilt-in-space chair with -5° to 40° tilt.




Ocean Dual VIP Ergo

An advanced tilt-in-space chair with -5° to 40° tilt and additional backrest recline.

NEW

Improved ergonomic headrest holder

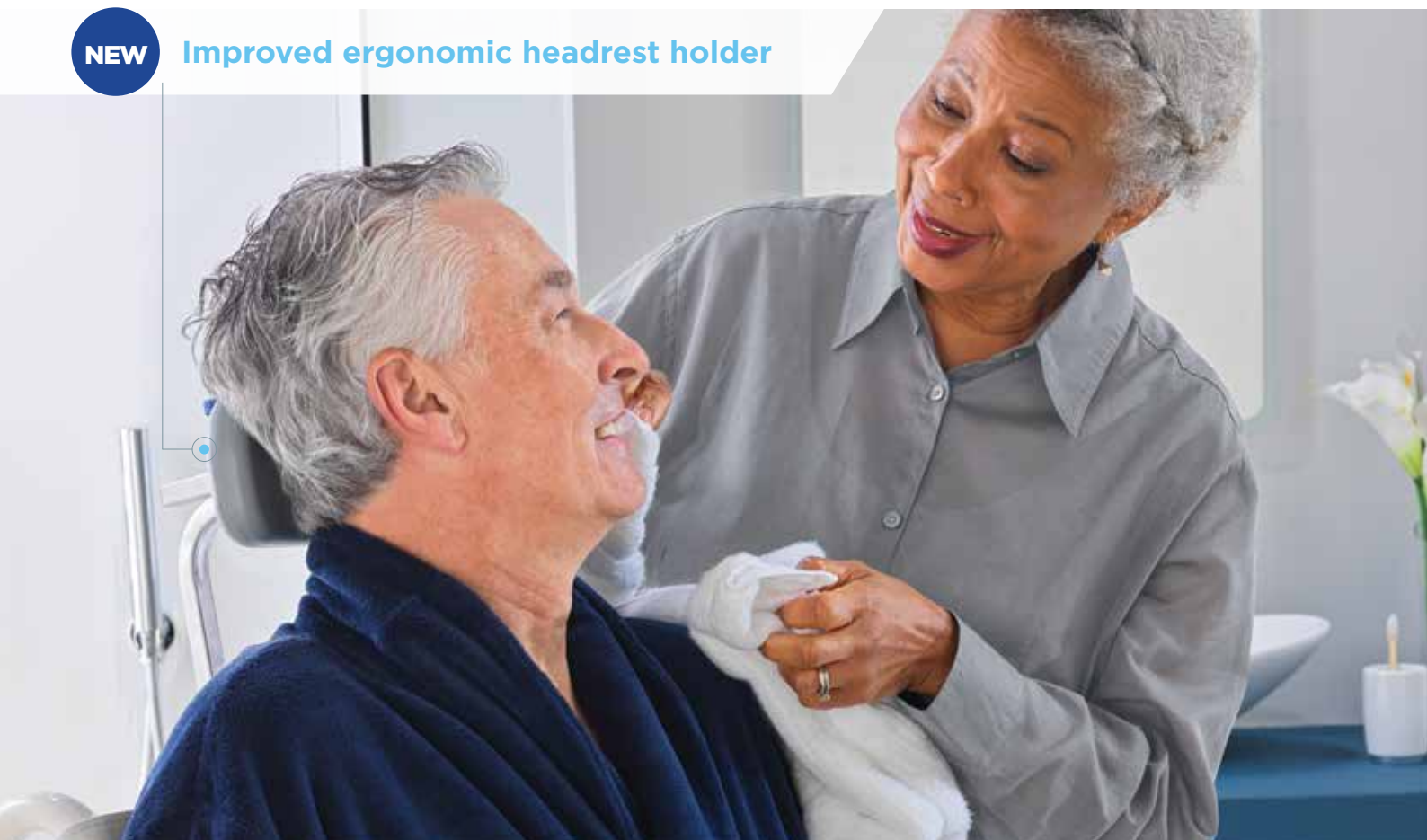

Our next generation of tilt-in-space shower chairs

The latest Aquatec Ocean VIP Ergo and Ocean Dual VIP Ergo are modular tilt-in-space shower chair commodes designed to suit a wide range of individuals and their postural needs. Both models benefit from a -5° to 40° seat tilt angle, whilst the Ocean Dual VIP Ergo has an additional backrest recline of 5° to 35° for those with more complex seating and positioning requirements.

As with all Ocean Ergo models, the VIP and Dual VIP Ergo can be tailored to each individual by making simple adjustments or by adding easy to fit, off the shelf accessories. The clinically-led design ensures easier use for the carer whilst providing optimum safety and comfort to the user.

NEW features

NEW

Improved ergonomic headrest holder

With multiple axis adjustment the ergonomic headrest holder offers **superior postural support for the head/neck region**. Adjustments are quick and easy by simply rotating the integrated hand-screws to achieve the correct position for the individual. Ideal for assessment and set up purposes with no tools required.

NEW

New tilt-stop function

An easy to apply optional function which **prevents the frame tilting anteriorly past 0°**. This ensures the frame remains horizontal to the ground and is an aid for those who are hoisted in/out of the chair.



Yes, you can.

Features and benefits



Unique tilt-in-space

The tilt-in-space mechanism means that during tilt, the chair moves downwards and forwards which keeps the user's centre of gravity close to that of the chair. This makes the tilt operation lighter and more intuitive for the carer and creates a smoother, more comfortable experience for the user. The tilt-in-space mechanism also means a reduced overall footprint without affecting its stability.



Contoured, ergonomic seat

Designed around the imprints people leave behind after being seated, to ensure a more upright and stable seating position for enhanced independence and comfort.

A 5° seat tilt lifts the knees into the natural squat position for toileting.

Toileting and transfer supports at the front to allow the user to position themselves into a comfortable position for toileting and showering.

Backrest recline

Backwards recline of between 5° - 35° is applicable for the Ocean Dual VIP Ergo.

5° - 35°

NEW

Improved armrest system

Newly designed armrests provide **increased support and comfort for the user**. The new and optional armrest pads offer individuals a softer surface area and resting position. Featuring an easy to operate and optional integrated safety lock system.

Adjustable armrests and footrests

Armrest height, footrest height and seat height can be easily adjusted.



Highly customisable

Available with an **extensive range** of positioning, toileting and seating accessories that can be used to support a wide range of conditions and postural needs.



Extremely stable

Enhanced stability for additional safety and reassurance, even at maximum tilt.



Smaller footprint

Optimised wheel base making it easier to manoeuvre in small spaces.



Steel frame

Stainless steel corrosion free frame.



Machine washable

Machine washable tension adjustable backrest.



Yes, you can.®

Tech Specs

For more comprehensive pre-sales information about this product, including the product’s user manual, please see your local Invacare website.



Seat width ▼



Seat depth ▼



Total width ▼



Total depth ▼



Width between armrests ▼

Ocean VIP Ergo	480 mm	480 mm	565 mm	1000 mm	475 mm
Ocean Dual VIP Ergo	480 mm	480 mm	565 mm	1000 mm	475 mm



Tilt angle seat ▼



Tilt angle back ▼



Weight ▼



Load capacity ▼

Ocean VIP Ergo	-5° - 40°	-	23.5 kg	150 kg
Ocean Dual VIP Ergo	-5° - 40°	5° - 35°	28.3 kg	150 kg



Frame colour



Yes, you can.®

rehab**hire** & sales

310 Lorimer Street Port Melbourne Victoria 3207

t: 03 9646 7200 | **e:** contact@rehabhire.com.au | www.rehabhire.com.au



FOR FURTHER INFORMATION:

New Zealand: Freephone 0800 468 222, email sales@invacare.co.nz or visit www.invacare.co.nz

Australia: Freephone 1800 460 460, email orders@invacare.com.au or visit www.invacare.com.au

Asia: Phone +61 2 8839 5330, email anzexportsales@invacare.com or visit www.invacare.com

Thailand: Phone +66 0 2821 5515 or visit www.invacare.com



Yes, you can.®