



Barwon Lift Recliner



Barwon Lift Recliner

The Barwon is a Lift-Recline Chair, the first choice for value in a high-quality Australian made lift chair. Choose from four standard chair size and an extensive range of feature adjustments to achieve a completely personalised fit and ensure comfortable usage. Every chair is fitted with carefully tested and selected electrical components and can be upholstered in your choice of premium fabric from Oscar's very own Versotela range.



FEATURES

- 130kg weight capacity
- Single motor wall saver
- Controller right side
- 9 volt battery
- 50mm casters
- Comes standard with Velcro attached back cushions
- Emergency battery backup
- Made to measure
- Available in wide range of Versotela fabrics
- 2 Year Warranty on mechanism, electrics and fabric
- 5 Year Warranty on timber frame
- 100% manufactured in Australia by Oscar Furniture

SPECIFICATIONS



	Seat Height	Seat Width	Seat Depth	Arm Height	Backrest Height	Scapular Height	Footrest Length
MINI SIZE							
MIN ADJUSTMENTS	-30mm	-150mm	-150mm	-60mm	-180mm	Linked to back height	-50mm
MEASUREMENT	456mm	540mm	455mm	199mm	636mm	429mm	N/A
MAX ADJUSTMENTS	+20mm	+50mm	+75mm	+60mm	+150mm	Linked to back height	+50mm
SIZE A							
MIN ADJUSTMENTS	-30mm	-150mm	-150mm	-60mm	-180mm	Linked to back height	-50mm
MEASUREMENT	476mm	540mm	490mm	180mm	647mm	440mm	N/A
MAX ADJUSTMENTS	+20mm	+50mm	+75mm	+60mm	+150mm	Linked to back height	+50mm
SIZE B							
MIN ADJUSTMENTS	-30mm	-150mm	-150mm	-60mm	-180mm	Linked to back height	-50mm
MEASUREMENT	496mm	540mm	520mm	190mm	688mm	481mm	N/A
MAX ADJUSTMENTS	+25mm	+50mm	+50mm	+60mm	+150mm	Linked to back height	+50mm
SIZE C							
MIN ADJUSTMENTS	-30mm	-150mm	-150mm	-60mm	-180mm	Linked to back height	-100mm
MEASUREMENT	546mm	540mm	545mm	190mm	718mm	511mm	N/A
MAX ADJUSTMENTS	+50mm	+50mm	+25mm	+60mm	+120mm	Linked to back height	+50mm

MODIFICATIONS

MEMORY FOAM

180KG WEIGHT CAPACITY

DUAL ACTION

SEAT WIDTH

SEAT HEIGHT

BACK HEIGHT

SEAT DEPTH

HEAD AND ARM REST COVERS

CONTROLLER LEFT SIDE

GLUTERAL RECESS

FOOTREST LENGTH

24 VOLT BATTERY

CUSHION CUTOUT

NO CASTERS

NO SHOWWOOD

INCREASE SIZE

TABLET ARM

ARM HEIGHT

MANUAL MECHANISM

LEG REST SIDE WALLS

GLUTERAL RECESS

REMOVABLE HEADREST



FABRIC SELECTION




VersoTela
Fabrics

All VersoTela fabrics come standard with a suite of high-end specifications, such as, antimicrobial coating, mold resistance, scratch resistance and fire resistance.

Our wide range of collections each boast unique properties that are carefully designed to suit every customer's needs. Whether it be built for use in a care facility, or as a functional but stylish addition to the living room, our VersoTela fabric options are sure to tick all the boxes when it comes to crafting your perfect chair.

Please get in touch with one of our friendly sales representatives for fabric samples and to find out more about which fabric is right for you.

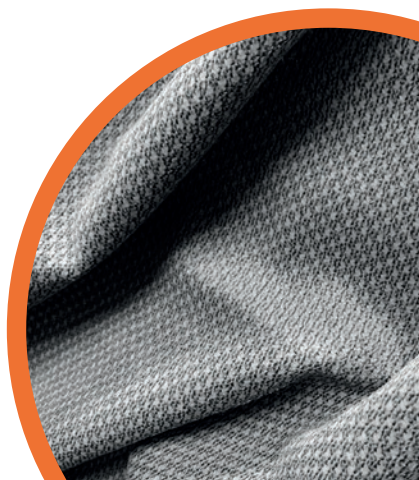



DESIGN


The Barwon comes standard with Velcro attached back cushions. The Velcro makes it easy to clean or replace cushions on site, and can help with minor adjustments to the cushion's position and fill.

Barwon three-roll back.

The top two rolls are filled with a fibrous filling that has a wonderful soft feel. Whilst the bottom roll is filled with a supportive foam. Barwon is the ultimate in softness and luxury.




ACTIONS

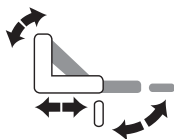

Single Motor (Wall-saver, Lift Assist)

The single motor wall saver operation of the Barwon allows the chair to move forward while the backrest reclines. This Oscar feature reduces the space that the chair occupies in the room because the reclined backrest does not hit the wall. The Lift, recline and footrest operate from one motor.

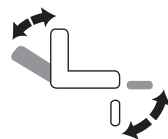

Dual Motor (Independent Recline)

Dual motor option lets the user recline the back without effecting the footrest, allowing a greater range of comfort.



Contact us today:
sales@oscarfurniture.com.au
03 5381 1404

6 King Drive,
Horsham VIC 3400

www.oscarfurniture.com.au

