

**Q-SERIES**

**THE  
POWER OF  
QUICKIE**



**AUSTRALIA**



Advanced engineering. Customised for you.



Q50 R Carbon



Q300



Q400



Q500/Q500 Mini



Q700



Q700-UP


# ADVANCED ENGINEERING. CUSTOMISED FOR YOU.

QUICKIE®, SEDEO®, JAY®, WHITMYER® and SWITCH-IT® bring you the most comprehensive range of premium, innovative and technically advanced powerchairs, each created individually for you whatever your clinical or lifestyle need.



# THE POWER OF QUICKIE

A complete power range from basic mobility to complex rehab.  
Whatever you need - without compromise.

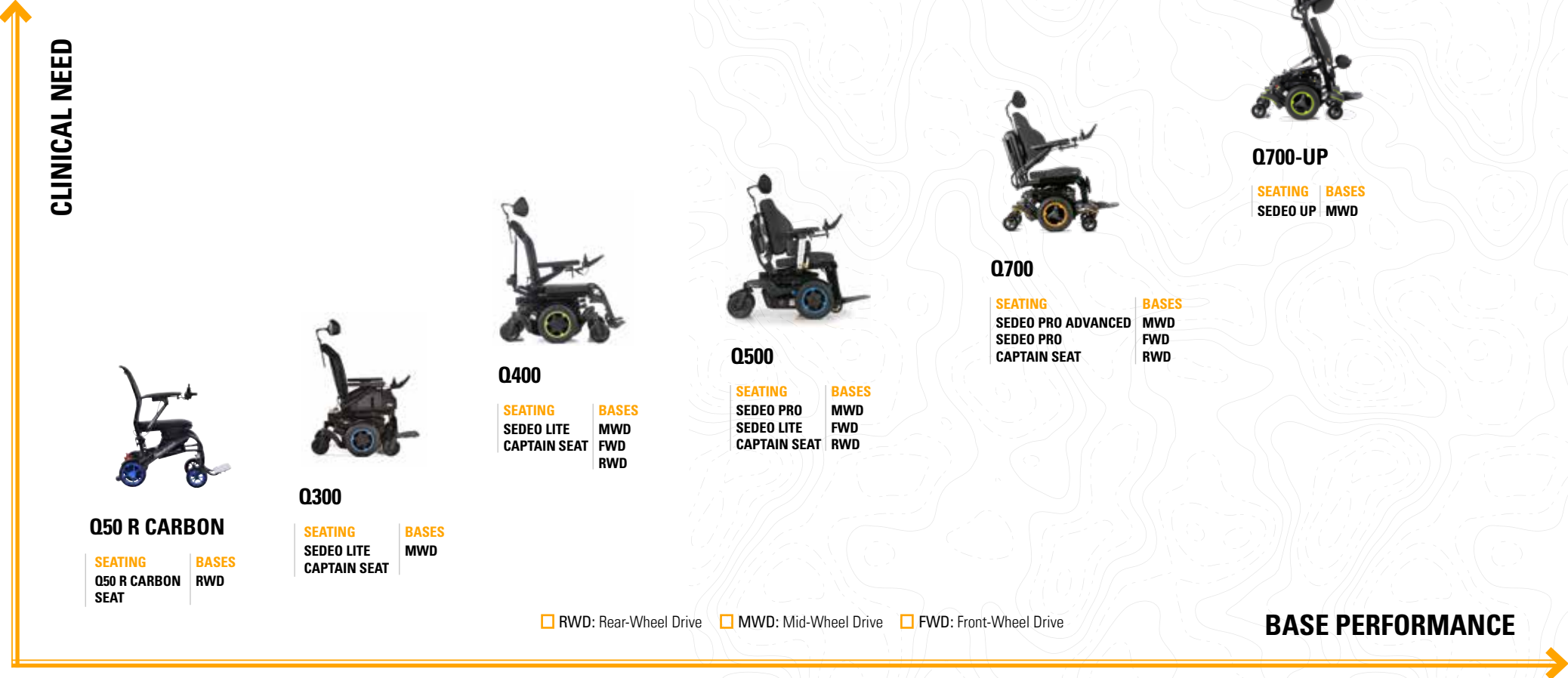

## SUMMARY

# THE NEW QUICKIE Q-SERIES POWER RANGE: A COMPLETE POWER PRODUCT PORTFOLIO UNDER ONE ROOF.

A wide range of bases, drive tractions, seating systems, electronics, controls and accessories, that can be combined to create the perfect powerchair for each individual need.



## SERIES

Q700

Q500

Q400

Q300

Q50 R CARBON

## BASES

Front-Wheel Drive

Mid-Wheel Drive

Rear-Wheel Drive

## SEATING SYSTEMS

SEDEO UP

SEDEO PRO ADVANCED

SEDEO PRO

SEDEO LITE

CAPTAIN SEAT

## CONTROLS

VR2

R-Net

SWITCH-IT Controls

Link-It Mounting System

Remote Seating App

## SUMMARY

# CHOICE WITHOUT COMPROMISE

The new QUICKIE Q-Series offers full flexibility to create the perfect powerchair for each individual user. Highly flexible configurations, interchangeable parts, quick and easy adjustments.



## BEST IN CLASS SEATING

The SEDEO seating range offers high-end seating systems to cover any clinical need. A choice of five seating systems: SEDEO UP, SEDEO PRO ADVANCED, SEDEO PRO and SEDEO LITE. All compatible with JAY cushions and backs and WHITMYER headrests to assure the best seating and positioning.



## BEST DRIVE BASES CHOICE IN THE MARKET

From cost effective compact bases to high performance outdoor bases, all with the most extensive choice of drive tractions. An impressive range of base levels, and drive tractions that can accommodate all sort of user preferences, lifestyle needs and pricing/reimbursement requirements.



## WITH MOST ADVANCED TECHNOLOGIES

Advanced engineering and innovative technologies applied to improve people's lives: Full Track Suspension, Memory Seating, advanced anti-shear backrests, biometrical repositioning, Remote Seating app (pressure management), SWITCH-IT controls and Link-It mounting systems and advanced R-Net Controls.



## EASIER TO LEARN

The common design and the modularity create a simplified user experience, that makes it very easy to learn, manage and prescribe.



## EASIER TO SERVICE

The modularity and flexibility of the bases, seatings and controls allows an easy adaptation to different users, making it perfect for refurbishments and allowing a great fleet management.

SERIES OVERVIEW

|                     | Q700-UP                                                                                             | Q700                                                                                                | Q500                                                                                                            | Q500 Mini                                                                                                       | Q400                                                                                                           | Q300                                                                                                         | Q50 R CARBON                                                                                          |
|---------------------|-----------------------------------------------------------------------------------------------------|-----------------------------------------------------------------------------------------------------|-----------------------------------------------------------------------------------------------------------------|-----------------------------------------------------------------------------------------------------------------|----------------------------------------------------------------------------------------------------------------|--------------------------------------------------------------------------------------------------------------|-------------------------------------------------------------------------------------------------------|
| Outdoor performance | <div><div></div><div></div><div></div><div></div><div></div></div>                                  | <div><div></div><div></div><div></div><div></div><div></div><div></div></div>                       | <div><div></div><div></div><div></div><div></div><div></div><div></div></div>                                   | <div><div></div><div></div><div></div><div></div><div></div><div></div></div>                                   | <div><div></div><div></div><div></div><div></div><div></div><div></div></div>                                  | <div><div></div><div></div><div></div><div></div><div></div><div></div></div>                                | <div><div></div><div></div><div></div><div></div><div></div><div></div></div>                         |
| Indoor performance  | <div><div></div><div></div><div></div><div></div><div></div><div></div></div>                       | <div><div></div><div></div><div></div><div></div><div></div><div></div></div>                       | <div><div></div><div></div><div></div><div></div><div></div><div></div></div>                                   | <div><div></div><div></div><div></div><div></div><div></div><div></div></div>                                   | <div><div></div><div></div><div></div><div></div><div></div><div></div></div>                                  | <div><div></div><div></div><div></div><div></div><div></div><div></div></div>                                | <div><div></div><div></div><div></div><div></div><div></div><div></div></div> Folds!                  |
| Suspension          | <div><div></div><div></div><div></div><div></div><div></div><div></div></div> <div>All Wheels</div> | <div><div></div><div></div><div></div><div></div><div></div><div></div></div> <div>All Wheels</div> | <div><div></div><div></div><div></div><div></div><div></div><div></div></div> <div>Drive wheels, Advanced</div> | <div><div></div><div></div><div></div><div></div><div></div><div></div></div> <div>Drive wheels, Advanced</div> | <div><div></div><div></div><div></div><div></div><div></div><div></div></div> <div>Drive wheels, Regular</div> | <div><div></div><div></div><div></div><div></div><div></div><div></div></div> <div>Drive wheels, Basic</div> | <div><div></div><div></div><div></div><div></div><div></div><div></div></div> <div>Drive wheels</div> |
| Motors              | 4 pole                                                                                              | 4 pole                                                                                              | 4 pole                                                                                                          | 4 pole                                                                                                          | 4 pole                                                                                                         | 4 pole                                                                                                       | 2 pole                                                                                                |
| Speed               | 8,10 kph                                                                                            | 8,10 kph                                                                                            | 10 kph                                                                                                          | 10 kph                                                                                                          | 10 kph                                                                                                         | 10 kph                                                                                                       | 6 kph                                                                                                 |
| Batteries           | 60 Ah / 72 Ah                                                                                       | 60 Ah / 72 Ah                                                                                       | 60 Ah                                                                                                           | 50 Ah                                                                                                           | 60 Ah                                                                                                          | 50 Ah                                                                                                        | 10 Ah Li-ion                                                                                          |
| Max. user weight    | 120 kg                                                                                              | 160 kg / 136 kg                                                                                     | 136 kg                                                                                                          | 136 kg                                                                                                          | 136 kg                                                                                                         | 136 kg                                                                                                       | 136 kg                                                                                                |
| Drive bases         | MWD / FWD                                                                                           | MWD / FWD / RWD /<br>PRO AVANCED                                                                    | MWD / FWD / RWD /                                                                                               | MWD                                                                                                             | MWD / FWD / RWD                                                                                                | MWD                                                                                                          | RWD / FOLDABLE                                                                                        |
| Seating             | SEDEO UP                                                                                            | SEDEO PRO ADVANCED/ SEDEO<br>PRO/ SEDEO LITE (136kg only) /<br>CAPTAIN SEAT (MWD 136kg only)        | SEDEO PRO / SEDEO LITE                                                                                          | SEDEO PRO /<br>CAPTAIN SEAT (MWD only)                                                                          | SEDEO LITE /<br>CAPTAIN SEAT (MWD only)                                                                        | SEDEO LITE /<br>CAPTAIN SEAT                                                                                 | COMFORT SEAT                                                                                          |
| Electronics         | R-Net                                                                                               | R-Net/VR2 CAPTAIN SEAT                                                                              | VR2/R-Net                                                                                                       | R-Net                                                                                                           | VR2                                                                                                            | VR2                                                                                                          | NVR2                                                                                                  |
| Crash Tested        | Yes                                                                                                 | Yes                                                                                                 | Yes                                                                                                             | Yes                                                                                                             | Yes                                                                                                            | Yes                                                                                                          | No                                                                                                    |

CREATED INDIVIDUALLY  
FOR YOUR CLINICAL AND  
LIFESTYLE NEEDS



## SERIES OVERVIEW

# QUICKIE Q-SERIES OVERVIEW

---

### Q700 SERIES: THE ULTIMATE EXPERIENCE

Chairs designed to offer a superb outdoor experience while allowing a good indoor use. A choice of three high performance drive bases with market leading suspension technologies, powerful motors and long lasting batteries to provide the ultimate outdoor performance. Combined with a choice of complex rehab seating systems (SEDEO UP STANDING WHEELCHAIR, SEDEO PRO ADVANCED & SEDEO PRO) to match users with specific or high clinical needs, and the most advanced electronics to access their world. Including also Captain Seat option.

### Q500 SERIES: MAKE YOUR OWN ADVENTURE

Robust outdoor powerchairs with a choice of three drive bases with advanced suspension technologies and powerful motors, for a superior outdoor experience. Combined with a choice of rehab seating systems (SEDEO PRO & SEDEO LITE) to provide great support and positioning to users from mild to medium clinical needs, and the most advanced electronics to access their world.

### Q500 MINI: ULTRA COMPACT PERFORMANCE SUITABLE FOR COMPLEX REHAB NEEDS

Ultra compact wheelchair with highly adjustable and customisable seating and electronics, a combination to match users with high clinical needs, but tight indoor environments.

### Q400 SERIES: ADVANCED OUTDOOR / INDOOR PERFORMANCE

Robust outdoor powerchairs with a choice of three drive bases, with regular motors and suspension for a great outdoor experience with advanced electronics. Available with SEDEO LITE seating systems and Captain Seat.

### Q300 SERIES: ULTRA COMPACT EVERY DAY PERFORMANCE

Entry level chairs with ultra-compact mid base for great indoor agility with the stability and traction for outdoor navigation and community access.

### Q50 R CARBON: TRAVEL LIGHTER

Ultra light weight foldable carbon frame for people on the go.  
14.5 kg and easily foldable for storage and transport.





## PERFORMANCE BASES

# DRIVE BASES OVERVIEW

The new QUICKIE powerchairs range offers the most complete choice of drive bases on the market, available in different series and with many seating systems.

The choice of a drive base affects how the chair performs in different situations and environments. The right choice for each user will depend on personal preferences, experience, clinical, lifestyle and environmental needs.

|                 | MWD       | FWD         | RWD         |
|-----------------|-----------|-------------|-------------|
| Manoeuvrability | ● ● ● ● ● | ● ● ● ● ● ● | ● ● ● ● ●   |
| Stability       | ● ● ● ● ● | ● ● ● ● ● ● | ● ● ● ● ● ● |
| Kerb climbing   | ● ● ● ● ● | ● ● ● ● ● ● | ● ● ● ● ● ● |
| Traction        | ● ● ● ● ● | ● ● ● ● ● ● | ● ● ● ● ● ● |
| Drive comfort   | ● ● ● ● ● | ● ● ● ● ● ● | ● ● ● ● ● ● |
| Ease of driving | ● ● ● ● ● | ● ● ● ● ● ● | ● ● ● ● ● ● |
|                 | MWD       | FWD         | RWD         |
| Q700 UP         | ✓         | ✓           | —           |
| Q700            | ✓         | ✓           | ✓           |
| Q500            | ✓         | ✓           | ✓           |
| Q400            | ✓         | ✓           | ✓           |
| Q300            | ✓         | —           | ✓           |

## PERFORMANCE BASES



### **MID-WHEEL DRIVE BASES (MWD)**

The mid-wheel drive bases have the best indoor performance. Their smallest turning circle provides great manoeuvrability. They are very intuitive and easy to drive, ideal for users with reduced hand control, driving with switches and users with poor spatial awareness.



### **FRONT- WHEEL DRIVE BASES (FWD)**

The front-wheel drive bases have great outdoor performance with good indoor manoeuvrability. Knee angle can be positioned at 90°. Maintains control on steep declines.



### **REAR-WHEEL DRIVE BASES (RWD)**


The rear-wheel drive bases are stable and have good kerb climbing ability. Intuitive drive that provides good control with non-proportional inputs.

# POWER SEAT FUNCTIONS




## C-ME® SEAT ELEVATION

No more waiting around or driving at reduced speed when elevated. C-Me enables you to elevate up to 12" while traveling up to 5 kmh. You can lift or descend at the touch of a button with no need to stop while transitioning. Just rise up and go at the same time!

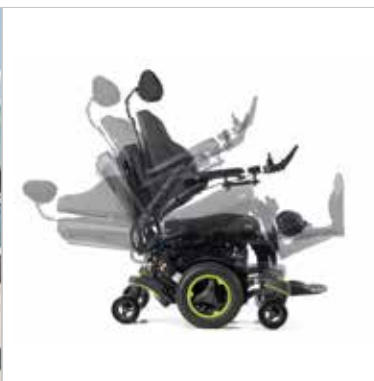

## POWER SHEAR REDUCTION RECLINE

Our power recline backrest offers 4" of shear reduction to minimize sliding and maintain good positioning.

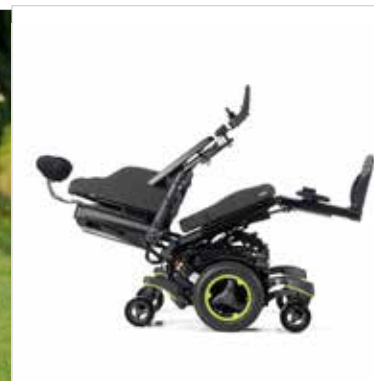

## POWER TILT

Center of gravity tilt providing pressure redistribution. You can stay in the chair for longer with more comfort and confidence knowing your skin is staying healthy.



## STANDING

Offering uncompromised standing position in a great performing mid-wheel-drive power wheelchair. Program the standing to meet your personal and functional needs.



## SEATING SYSTEMS

# SEATING OVERVIEW

|                                         | SEDEO UP                    | SEDEO PRO<br>ADVANCED     | SEDEO PRO                          | SEDEO LITE    | CAPTAIN SEAT   | Q50 R CARBON |
|-----------------------------------------|-----------------------------|---------------------------|------------------------------------|---------------|----------------|--------------|
| <b>User clinical need</b>               | High / Standing             | High                      | Medium to High                     | Medium to Low | Low            | Low          |
| <b>Seat width</b>                       | 400 - 500 mm                | 400 - 560 mm              | 380 - 560 mm                       | 385 - 510 mm  | 16 - 20 inches | 17 inches    |
| <b>Seat depth</b>                       | 400 - 540 mm                | 400 - 560 mm              | 400 - 560 mm                       | 410 - 560 mm  | 16 - 22 inches | 17 inches    |
| <b>Backrest height</b>                  | 560 - 610 mm                | 560 - 700 mm              | 440 - 620 mm                       | 460 - 510 mm  | 25 inches      | 17 inches    |
| <b>Configurability</b>                  | Low                         | Low                       | High                               | Medium        | —              | Foldable     |
| <b>Standing function</b>                | ✓                           | —                         | —                                  | —             | —              | —            |
| <b>Lay down position</b>                | ✓                           | ✓                         | —                                  | —             | —              | —            |
| <b>Programmable seating</b>             | ✓                           | ✓                         | —                                  | —             | —              | —            |
| <b>Powered recline</b>                  | 90° - 175°<br>biomechanical | 90°-175°<br>biomechanical | 90°-130°/90°-150°<br>biomechanical | 85°-120°      | —              | —            |
| <b>Manual recline</b>                   | —                           | —                         | 90°-120°                           | 90°-120°      | 90°-115°       | —            |
| <b>Anti-shear back movement</b>         | 180 mm                      | 180 mm                    | 100 mm (optional)                  | —             | —              | —            |
| <b>Pressure release reminders</b>       | ✓                           | ✓                         | —                                  | —             | —              | —            |
| <b>Tilt angle</b>                       | 22°                         | 50°                       | 50°                                | 30°/50° *     | 30°/50° *      | —            |
| <b>Anterior Tilt</b>                    | 22°                         | 25°                       | N/A                                | N/A           | N/A            | —            |
| <b>Lift</b>                             | 200 mm                      | 300 mm                    | 300 mm                             | 300 mm        | 300 mm         | —            |
| <b>JAY &amp; WHITMYER compatibility</b> | ✓                           | ✓                         | ✓                                  | ✓             | JAY Cushion    | JAY Cushion  |
| <b>Remote Seating app</b>               | ✓                           | ✓                         | —                                  | —             | —              | —            |
| <b>Max. user weight</b>                 | 120 kg                      | 160 kg                    | 160 kg                             | 136 kg        | 136 kg         | 136 kg       |
| <b>Crash Tested</b>                     | Yes, ISO                    | Yes, ISO                  | Yes, ISO                           | Yes, ISO      | Yes, ISO       | No, Foldable |

## SEATING SYSTEMS



### SEDEO UP

High-end biomechanical seating, with many programmable seating positions and standing function. All the benefits of SEDEO seating with additional standing function. SEDEO UP provides a standing angle up to 85° with a highly rigid and stiff standing module.



### SEDEO PRO ADVANCED

High-end biomechanical seating, with many programmable seating positions. SEDEO PRO ADVANCED covers the most demanding clinical needs with 180 mm anti-shear back movement, biometric repositioning function, infinitely adjustable seat and back, easy entry system, 90° knee angle in seated position and pressure relief reminders.




### SEDEO PRO

High-end fully adjustable and configurable seating frame with a wide range of comfort cushions and positioning options. Aimed at meeting moderate to complex clinical needs. With a wide range of armrests, armpads, legrests, seat and back cushions, headrests and positioning accessories, SEDEO PRO offers great flexibility to adapt to the individual user needs.




### SEDEO LITE

Multi-adjustable seat frame with the option of either SEDEO LITE backrest cushions, Jay backrests or sling back upholstery. SEDEO LITE is designed for users with medium to low clinical needs that require adjustable seating system. The backrest, armrests and legrests can be adjusted to fit the user. Includes seat frame with seat rail for additional mounting solutions.

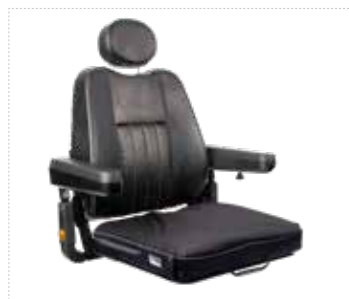

### CAPTAIN SEAT

Provides a relaxed sitting position and features a semi-recline, fold-down backrest, adjustable seat size, and convertible seat pan that allows you to use your preferred cushion. Available with power lift or 30°/50° tilt and power centermount legrest.

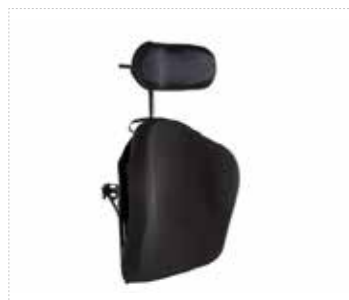

### WHITMYER & JAY SEATING

Compatible with all SEDEO seating systems, the JAY family of cushions and backs offers superior clinical seating for pressure management, positioning and comfort. WHITMYER headrests are also compatible with all SEDEO seating systems offering the best head support and positioning.



### Q50 R CARBON

Simple comfortable seating for people on the go. Backrest has contoured and breathable fabric. Seat cushion is designed for stable support.

## CONTROLS & CONNECTIVITY

# VR2 & R-NET

## ADVANCED JOYSTICKS

The QUICKIE Q-Series wheelchairs come with the reliable VR2 control or the superior R-Net advanced joystick. A large, easy-to-read HD display and simple to control paddle switches allow for intuitive operation. (Depends on model).

### RELIABLE

Industry standard VR2 and R-Net systems offer easy to use and highly durable drive controls.

### PROGRAMMING

VR2 and R-Net both provide advanced programming to customise the chair to your requirements and environments.

### CONNECT

Seamlessly take control of your environment using the integrated Infrared functionality. Control your home entertainment system, TVs, DVDs, docking stations and other home entertainment systems simply through your joystick (R-Net only).

### BLUETOOTH & IDRIVE

Control your phone, computer, tablet and any other Bluetooth and i-devices with this useful optional extra. Connect and control up to four devices (R-Net only).

### QUICKIE'S PATENTED ASSIGNABLE BUTTONS

Change the function on any button on your joystick to any command you desire, thanks to QUICKIE's patented assignable buttons. You can also assign a second (long press) function to buttons, giving you quicker, simpler and smoother access to functions you commonly use without having to use your joystick (R-Net only).



**R-Net**  
Q500, Q500 Mini,  
Q700, Q700 UP




**VR2**  
Q400, Q300

## ALTERNATIVE CONTROLS

# SWITCH-IT INNOVATIVE CONTROLS

---

Conveniently access your 5 programmable positions with the Ctrl+5 – a perfect combination with Memory Seating (Memory Seating with SEDEO PRO ADVANCED and SEDEO UP only). It allows you to go straight to the function you need at the touch of a button (R-Net only).

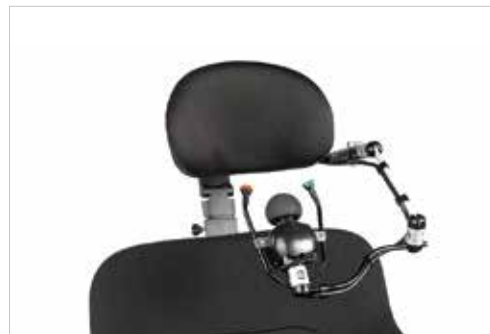

If driving with a traditional joystick is not the best option, alternative drive controls allow you to fully experience your chair's driving capabilities and benefits.

# LINK-IT MOUNTING SYSTEM

---

Need to use additional switches or special controls? No problem, easily connect them using the universal Link-It mounting system.





## DESIGN AND APPEARANCE

# DISTINGUISHED STYLE

With the Q-Series, we've taken style to the next level. With a matte-black chassis finish, hidden wiring and components combined with sharp-tailored upholstery, it's a sleek and a contemporary design with a variety of colours to choose from.

### FANCY A CHANGE?

A choice of contemporary colour inlays to perfectly complement its striking matt-black bodywork.

Simply replace the colour inlays. Styling your powered wheelchair has never been so simple.



\* The number of colours available may differ depending on the seat and base combination.





Sunrise Medical Australia Pty Ltd.  
11 Daniel Street  
Wetherill Park NSW 2164  
Phone: (02) 9678 6600  
Email: [enquiries@sunrisemedical.com.au](mailto:enquiries@sunrisemedical.com.au)  
[www.SunriseMedical.com.au](http://www.SunriseMedical.com.au)

© 2023 Q-Series Range Brochure AU Rev E

For further information on the full specification, options and accessories, please refer to the order form. All information is subject to change without notification. Please consult **Sunrise Medical** with any queries you may have.

Pictures and illustrations may be overseas models and some seating and accessories not available in Australia.

If you are visually impaired, this document can be viewed in PDF format at [SunriseMedical.com.au](http://SunriseMedical.com.au).

