

rehabhire & sales

YOUR  
LEADING  
HEALTHCARE  
EQUIPMENT  
PARTNER

novis

# AreaCare® Range

A comprehensive range of static foam pressure relieving surfaces designed to maintain skin integrity and help prevent pressure injury, while delivering maximum comfort.



1300 000 030 | [rehabhire.com.au](http://rehabhire.com.au)


YOUR  
LEADING  
HEALTHCARE  
EQUIPMENT  
PARTNER

## AreaCare Range

As specialists in the supply of pressure care and manual transfer equipment solutions designed to help prevent pressure wound and skin integrity issues, Novis support a range of clinical healthcare products designed for hospitals, aged care facilities, rehabilitation providers and the home care community.

With a long history delivering effective pressure care management solutions, we are committed to delivering successful patient outcomes through product innovation, professional collaboration and clinical excellence.


Our product portfolio has been engineered with the goal of enhancing the comfort and safety of patients; reducing the time and concerns of care staff; and improving clinical outcomes for care facilities.

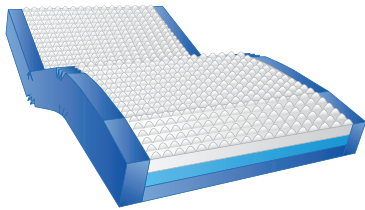

### AreaCare Optimum

- Castellated, multi-density, three layer, zoned design for maximum comfort
- Immersive soft memory foam top layer
- Hinged side walls which aid bed profiling and help prevent falls
- Heel Slope zone for localised pressure care


	Single	King Single
Capacity	220kg	220kg

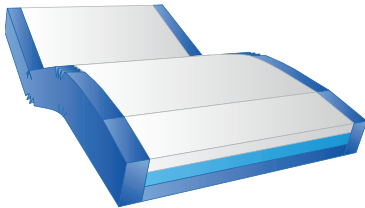

### AreaCare Premium

- Multi-density, three layer, zoned design for comfort and support
- Supportive firm memory foam top layer
- Hinged side walls which aid bed profiling and help prevent falls
- Heel Slope zone for localised pressure care


	Single	King Single
Capacity	220kg	220kg

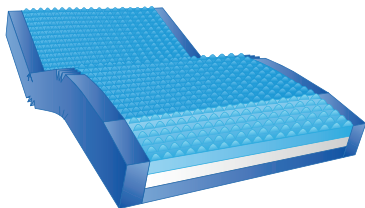

### AreaCare Comfort

- Castellated, high resilience non memory foam top layer for long term use
- Memory foam heel section for maximum focused immersion
- Hinged side walls which aid bed profiling and help prevent falls
- Heel slope zone for localised pressure care

	Single	King Single
Capacity	190kg	190kg

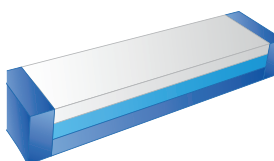

### AreaCare Extension

- 200mm extension for beds with sleeping platform extension
- Memory foam top layer provides additional pressure relief under 'at risk' heel zone

	Single	King Single
Capacity	NA*	NA*

\*Same capacity as mattress used

## An updated range with extensive features to offer every patient optimum care and comfort

Multi-density foam layers for maximum strength and comfort – made from premium material, each foam layer is made of varying densities to maximise the redistribution of pressure.

### Top Cover

Latex free, multi-stretch polyurethane coated top cover for reduced shear and friction. Vapour permeable, bacteriostatically treated, fire retardant and water resistant. Double-layered waterfall skirt provides protection from moisture ingress.

Cover includes an asset management panel and all relevant TGA required information.

### Top Layer

#### AreaCare Optimum

Features a castellated memory foam top designed to enhance pressure redistribution and articulation with movement.

#### AreaCare Premium

Supportive firm memory foam top layer, this foam contours to the body's shape, promoting even pressure distribution.

#### AreaCare Comfort

Features a castellated non memory foam top layer with higher resilience, allowing for easier repositioning.

### Infection Control

Top cover seams are welded and/or tape sealed for maximum infection control.

All models feature a 360° double-layered waterfall skirt to direct fluids away from mattress base.

### Inner Cotton Cover

Foam covered by a removable inner cotton cover for added foam protection and reduced shear.

### Support Layer

High resilience foam for maximum pressure redistribution.

### Hinged Base

Six hinge system to conform to profiling beds and enhance mattress flexibility.

### Core Layer


Firm core layer and edges provide form, strength and stability for maximum safety and reduce risk of falls.

### Side Bolster

All models feature a side bolster to provide a firm surface for reduced risk of falls when sitting on edge of mattress.

### Base Cover

Hard-wearing, welded and water resistant base cover material with 360 zipper for easy removal and two carry handles per side.

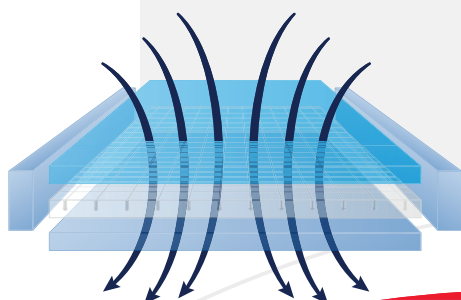

AreaCare Climate offers a unique airflow design to improve microclimate.





	AreaCare Optimum   Soft		AreaCare Premium   Firm		AreaCare Comfort	
Size	Single	King Single	Single	King Single	Single	King Single
Length	1980mm	1980mm	1980mm	1980mm	1980mm	1980mm
Width	880mm	1050mm	880mm	1050mm	880mm	1050mm
Height	150mm	150mm	150mm	150mm	150mm	150mm
Safe Working Load	220kg	220kg	220kg	220kg	190kg	190kg
Layers	3					
Side Walls	✓					
360° Zipper	✓					
Hinged	✓					
Waterfall Skirt	✓					
Heel Slop Taper Zone	✓					
Warranty	Body: 6 years   Cover: 2 years					
Castellation   Convolution	✓	✓	✗	✗	✓	✓
High Resilience PU Foam	✗	✗	✗	✗	✓	✓
Code	BSPS6003-S	BSPS6003-KS	BSPS6002-S	BSPS6002-KS	BSPS6001-S	BSPS6001-KS

Notes

Lined area for notes.

# rehabhire&sales

1 300 000 030

320 Lorimer Street Port Melbourne  
contact@rehabhire.com.au  
rehabhire.com.au

# rehabhire&sales **REPAIRS AND SERVICE**

REGULAR MAINTENANCE OR REPAIRS

1 300 000 030

Unit 6 | 52 Wirraway Drive Port Melbourne  
service@rehabhire.com.au  
rehabhire.com.au/repairs-and-service

# rehabinstallation

CAN MODIFY EVERY ROOM IN YOUR HOME TO  
ENSURE IT IS ACCESSIBLE, SAFE, AND COMFORTABLE

1 300 000 030

320 Lorimer Street Port Melbourne  
contact@rehabinstallation.com.au  
rehabinstallation.com.au



Registered NDIS Provider