

rehabhire & sales

YOUR
LEADING
HEALTHCARE
EQUIPMENT
PARTNER

A&D
A&D Weighing

Multi-Function Professional Blood Pressure Monitor

UM-211



1300 000 030 | rehabhire.com.au

YOUR
LEADING
HEALTHCARE
EQUIPMENT
PARTNER

The A&D Medical Professional Blood Pressure Monitor is clinical grade and simple to use, designed to deliver clinicians the highest level of precision and accuracy required to measure blood pressure. This economical blood pressure monitor is perfect for pharmacies, nursing homes, clinics, hospitals, other medical office environments, and home care.

Features and Benefits

- Dual measurement mode (Auscultatory and Oscillometric)
- Cuff holder and grip for carrying
- Durable and chemical resistant body
- Five optional cuff sizes available
- Large LCD with backlight
- Rechargeable battery
- Adjustable inflation pressure (AUTO/220/250/280)
- Irregular HeartBeat (IHB) indicator
- Room temperature display



Stand-by display
(Room temperature)



Measurement display
(with backlight)



In combination
with optional
mobile stand


What is in the Box

- 1 Blood Pressure Monitor
- 1 Blood Pressure Cuff
- 1 Nickel-Metal Hydride Battery
- 1 AC Adapter

UM-211's auscultatory mode is designed for daily and frequent use. Selectable deflation speed (2.5mmHg / sec. or 5.0mmHg / sec.) and quick start for auscultatory mode are features unique to A&D.



Outpatient

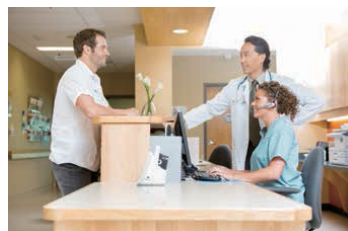

Stored and charged
at nurse station



Easy to carry



Measurement at ward

Technical Specifications

Measurement Method	Oscillometric and Auscultatory
Measurement Range	Pressure: 0 to 299mmHg Systolic Pressure: 60 to 279mmHg Diastolic Pressure: 40 to 200mmHg Pulse: 40 to 200 beats/minute
Measurement Accuracy	Pressure: ± 3 mmHg Pulse: $\pm 5\%$
Temperature Accuracy	$\pm 2.5^{\circ}\text{C}$ (36.5°F) ($+5^{\circ}\text{C}$ to $+40^{\circ}\text{C}$)
Temperature Unit	$^{\circ}\text{C}$ or $^{\circ}\text{F}$
Battery Life	Approximately 300 Measurements
Clinical Test	According to ISO81060-2 2013 IEC 60601-1-2:2007
Memory	99 Measurements
Power Supply	Oscillometric and Auscultatory
Dimensions (WxHxD)	Approximately W: 120mm H: 200mm D: 140mm
Unit Weight	Approximately 650g excluding batteries

Optional Accessories

Product	Code	Measurement Arm Circumference
Extra Small Cuff	UM-CUF-KS-SS	120mm - 170mm
Small Cuff	UM-CUF-KS-SA	160mm - 240mm
Adult Cuff	UM-CUF-KS-A	220mm - 320mm
Large Adult Cuff	UM-CUF-KS-LA	310mm - 450mm
Extra Large Cuff	UM-CUF-KS-LL	410mm - 500mm
Rolling Stand with Basket	UM-211-ST UM-21x-ST-OR27	



Notes

rehabhire&sales

1 300 000 030

320 Lorimer Street Port Melbourne
contact@rehabhire.com.au
rehabhire.com.au

rehabhire&sales **REPAIRS AND SERVICE**

REGULAR MAINTENANCE OR REPAIRS

1 300 000 030

Unit 6 | 52 Wirraway Drive Port Melbourne
service@rehabhire.com.au
rehabhire.com.au/repairs-and-service

rehabinstallation

CAN MODIFY EVERY ROOM IN YOUR HOME TO
ENSURE IT IS ACCESSIBLE, SAFE, AND COMFORTABLE

1 300 000 030

320 Lorimer Street Port Melbourne
contact@rehabinstallation.com.au
rehabinstallation.com.au



Registered NDIS Provider