

rehabhire&sales

YOUR
LEADING
HEALTHCARE
EQUIPMENT
PARTNER

A&D
A&D Weighing

Digital Blood Pressure Monitor


UA-767S-W



1300 000 030 | rehabhire.com.au

YOUR
LEADING
HEALTHCARE
EQUIPMENT
PARTNER

The UA-767S-W (which includes a wide range cuff) sets the standard for personal use. This monitor is easy to use and features the %IHB function with memory of 60 readings.

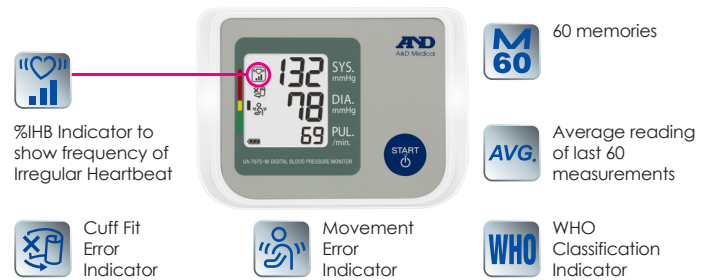

About %IHB

Features

The frequency of IHB (Irregular Heart Beat) detection is divided into different grades and shown by icons.



0-24%: No Indication
25-49%: Grade 1
50-74%: Grade 2
75-100%: Grade 3

Benefits

You can see how often IHB has been detected (indicating the level of risk). If the grade is high, you should consult a doctor immediately. This %IHB function is not to be used for diagnosis, but for monitoring purposes only.

Specifications

Measurement Method	Oscillometric
Measurement Range	Pressure: 0 to 299mmHg Systolic Pressure: 60 to 279mmHg Diastolic Pressure: 40 to 200mmHg Pulse: 40 to 180 beats/minute
Measurement Accuracy	Pressure: ± 3 mmHg Pulse: $\pm 5\%$
Display Type	Digital LCD Pressure/pulse displayed simultaneously
Number of Measurements	Approx. 700 times LR6 (alkaline batteries) Approx. 200 times R6P (manganese batteries) With pressure value 180 mmHg, room temperature 23 °C.
Memory Storage	60
Memory Data Average	Maximum 60 data average of systolic, diastolic, and pulse values
Operating Conditions	+10°C to +40°C 15 %RH to 85 %RH 800 to 1060hPa
Transport Storage Conditions	-20°C to +60°C 10 %RH to 95 %RH
Power Supply	4 x 1.5V batteries (R6P, LR6 or AA) or AC adapter (ECTB-182-O: not included)
Dimensions	Approx. W: 140mm H: 60mm D: 105mm
Unit Weight	Approx. 245g excluding batteries

Optional Accessories

- Large SlimFit Cuff | 310mm-450mm: UA-CUFBOXLA-C
- AC Adaptor: ECTB-182-O

rehabhire&sales

1 300 000 030

320 Lorimer Street Port Melbourne
contact@rehabhire.com.au
rehabhire.com.au

rehabhire&sales **REPAIRS AND SERVICE**

REGULAR MAINTENANCE OR REPAIRS

1 300 000 030

Unit 6 | 52 Wirraway Drive Port Melbourne
service@rehabhire.com.au
rehabhire.com.au/repairs-and-service

rehabinstallation

CAN MODIFY EVERY ROOM IN YOUR HOME TO
ENSURE IT IS ACCESSIBLE, SAFE, AND COMFORTABLE

1 300 000 030

320 Lorimer Street Port Melbourne
contact@rehabinstallation.com.au
rehabinstallation.com.au



Registered NDIS Provider