

rehabhire & sales

YOUR
LEADING
HEALTHCARE
EQUIPMENT
PARTNER

novis

ProCair Bariatric

Designed for hospital and aged care clients with higher weight needs and very high risk of pressure injuries.



1300 000 030 | rehabhire.com.au

YOUR
LEADING
HEALTHCARE
EQUIPMENT
PARTNER

ProCair Bariatric

The ProCair Bariatric Alternating Pressure Mattresses are a practical solution for bariatric patients at risk of pressure injuries delivering effective care without compromising on comfort or affordability. With a 450kg Safe Working Load, advanced 1-in-4 air cell alternation, and a fully sealed base for enhanced infection control, ProCair Bariatric provides reliable, high-performance pressure care for hospital and community environments.



Safe Working Load

Designed to safely support patients up to 450kg for effective bariatric care.

Top Cover

Industry leading Carrflex Extreme fabric reduces shear and friction forces while maintaining the highest levels of water resistance and breathability while maintaining the ideal microclimate.

Bed Straps

Designed to safely support patients up to 450kg for effective bariatric care.

Patient Records Panel

Waterfall Skirt

The waterproof top cover has a full length waterfall skirt and welded seams to provide maximum control against infection.

Umbilical Cap

For 24 hour transport mode after 3 hour active therapy from backup battery.

AutoCair


No more guesswork in pressure settings. AutoCair's automatic weight-sensing technology continuously adjusts pressure levels to keep the client within a safe therapeutic range—aligned with the International Clinical Practice Guideline (2025).

Fowler Boost

Auto Fowler Boost mode adjusts mattress pressure to provide increased support under the sacral area when the backrest is raised preventing risk of 'bottoming out' under the increased load.

Back Up Battery

An inbuilt Back Up Battery unit provides the re-assurance of emergency back up power in event of a power failure or accidental power disconnection. The pump will continue to operate using reserve battery power for continuous operation of alternation of up to 3 hours. Also ideal for alternating therapy during extended patient transportation periods.

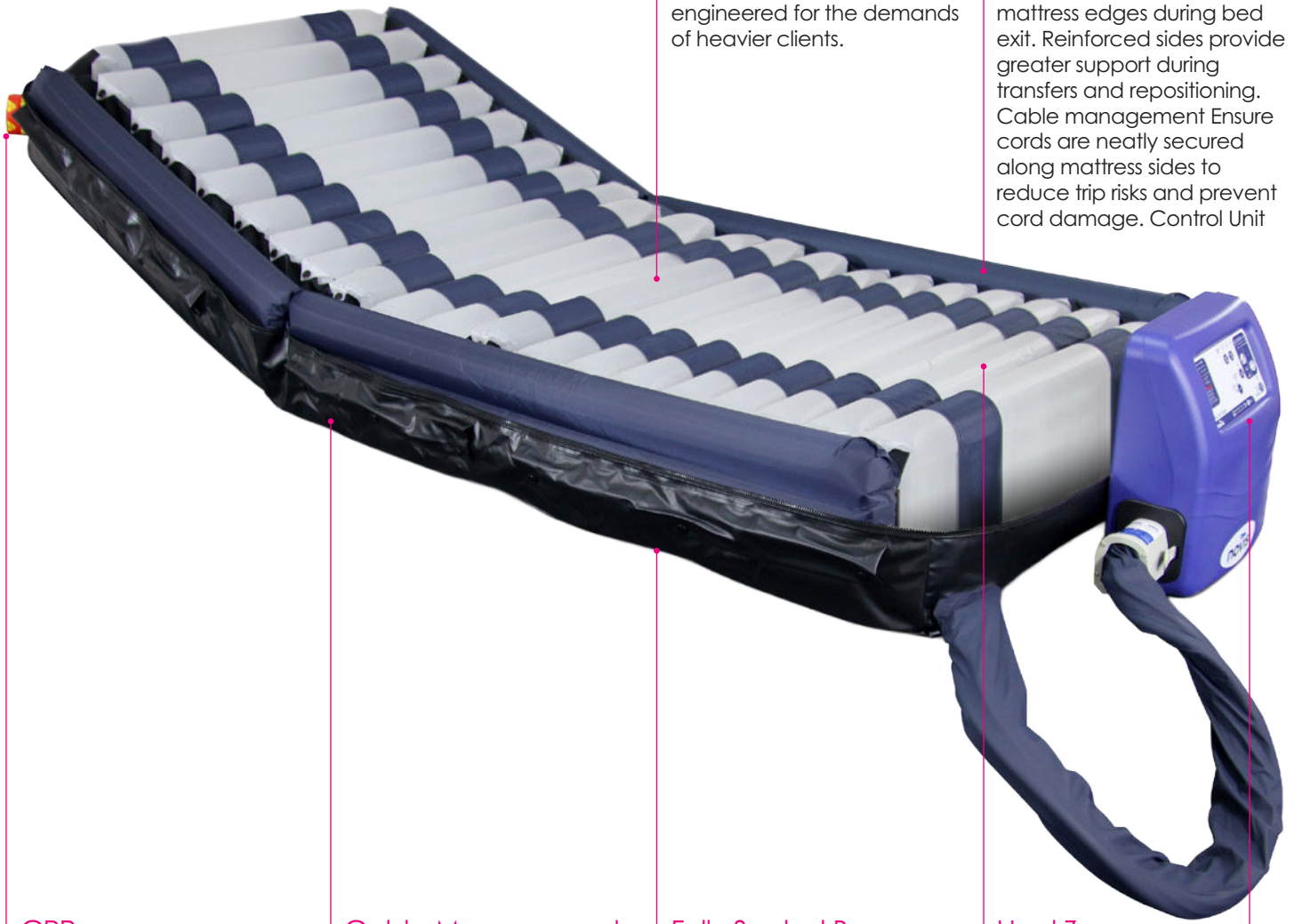

Cells

Lightweight and flexible pure TPU cells allow for maximum immersion. Deep cell design is engineered for the demands of heavier clients.

Static Side Bolsters

Enhance patient safety and stability by minimising the risk of cell collapse at the mattress edges during bed exit. Reinforced sides provide greater support during transfers and repositioning. Cable management Ensure cords are neatly secured along mattress sides to reduce trip risks and prevent cord damage. Control Unit

CPR

Quick-release CPR access at the head of the mattress allows rapid deflation in emergencies.

Cable Management

Ensure cords are neatly secured along mattress sides to reduce trip risks and prevent cord damage.

Fully Sealed Base

For maximum infection control.

Heel Zone

Individually detachable heel zone cells to allow complete off-loading of the calcaneus.



Control Unit

Dual compressor output allows for superior 450kg Safe Working Load. Equipped with a back up battery and automatic weight sensing technology.

Specifications

Size	King Single	Wide
SKU	APMPC-B04S	APMPC-B04SW
Length	2000mm	2000mm
Width	1050mm	1219mm
Height	250mm	250mm
Weight	16.5kg	18.5kg
Risk Category	High to Very High	
Fowler Boost	Automatic	
SWL	40kg - 450kg	
Total Number of Cells	20	
Static Head Cells	3	
Alternation Cycle	1:4	
Cell Type	Deep Cell	
Cell Material	Pure TPU	
Top Cover	Carrflex Extreme	
Control Unit	External with built-in hooks	
Compliance	ARTG 289458	
Warranty	3 years	

rehabhire & sales

1 300 000 030

320 Lorimer Street Port Melbourne
contact@rehabhire.com.au
rehabhire.com.au

rehabhire & sales **REPAIRS AND SERVICE**

REGULAR MAINTENANCE OR REPAIRS

1 300 000 030

Unit 6 | 52 Wirraway Drive Port Melbourne
service@rehabhire.com.au
rehabhire.com.au/repairs-and-service

rehabinstallation

CAN MODIFY EVERY ROOM IN YOUR HOME TO
ENSURE IT IS ACCESSIBLE, SAFE, AND COMFORTABLE

1 300 000 030

320 Lorimer Street Port Melbourne
contact@rehabinstallation.com.au
rehabinstallation.com.au



Registered NDIS Provider