

rehabhire & sales

YOUR
LEADING
HEALTHCARE
EQUIPMENT
PARTNER



Blake
Medical Distribution

Geo-Matrix™ Vinga




1300 000 030 | rehabhire.com.au

The **Geo-Matrix™ Vinga** cushion perfectly combines the pressure and shear reduction of Geo-Matrix™ gel with the plush comfort of gel infused memory foam.

The Vinga was designed to maximise comfort and protection while maintaining a stable posture. With its multi foam core and the firm ischial shelf the Vinga finds the balance between comfort and stability.


The Vinga is designed to be light-weight and maximise comfort, while maintaining stable posture.

Each Vinga cushion comes complete with either a breathable or soft touch incontinent cover, as well as a 4 way stretch protective inner cover.

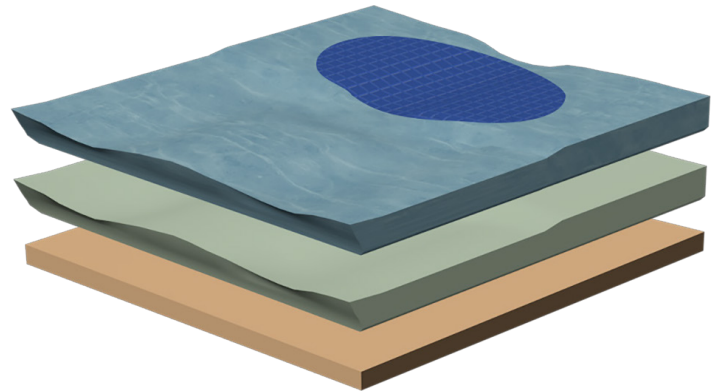

Features and Benefits

- Cushion remains 2-3 degrees cooler under seated body temperature to assist with microclimate management
- Minimises shear through column segmentation
- Geo-Matrix™ Gel does not change structure in climates of extreme heat nor extreme cold

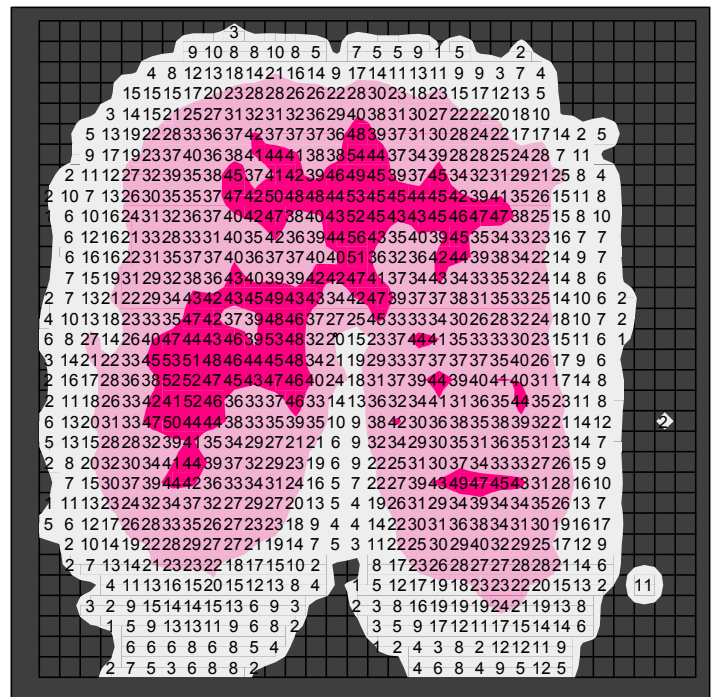

Client: Female 53kg

Specifications

Construction

Closed cell foam base and lateral stabilising wedges. High resiliency foam for the base, and ischial shelf. Core is comprised of a contoured injected foam. The top layer features a gel infused memory foam at the front of the cushion and a zoned Pillow Top Geo-Matrix™ gel insert along the back half of the cushion.

All of Blake Medical products are free from phthalates and latex. This includes all medical fabrics, Geo-Matrix™ gel, moulded polyurethane, standard polyurethane and specialty foam.

Covers

Two outer cover options: Breathable Spacetex fabric with incontinence protection OR Elastatex soft touch outer incontinence cover.

Inner cover: Elastatex breathable, wipe-able, antibacterial, fluid proof cover.

Compliance

All cushions meet EN-1021 part 1 & 2 requirements for fire retardancy.

Weight Capacity

Recommended Maximum 125kg

Warranty

24months

Advantages of Geo-Matrix™ Gel



Durable



Comfortable



Cleanable



Friction & Shear Reducing



Increases Sit Times



Microclimate Management



Pressure Redistributing



Vibration Dampening



Won't Puncture



Washing Instructions

Covers: 60°C in washing machine. Hang to air dry. Do not bleach.

Geo-Matrix™ gel and foam: Wipe down gel using any cleaner that does not contain bleach and allow to air dry.



1300 000 030

320 Lorimer Street Port Melbourne
contact@rehabhire.com.au
rehabhire.com.au



1300 000 030

Unit 6 | 52 Wirraway Drive Port Melbourne
service@rehabhire.com.au
rehabhire.com.au/repairs-and-service



1300 000 030

320 Lorimer Street Port Melbourne
contact@rehabinstallation.com.au
rehabinstallation.com.au